

home care bed Arminia III



home care bed Arminia III

> performance

4-section mattress base 90 x 200 cm, mattress surface height adjustment of approximately 40 – 80 cm, back rest and thigh rest as well as the mattress surface height are adjustable by electric motors. The foot section height can be lowered by tilting the mattress base. Metal parts are sealed with a polyester powder lacquer coating, in the colour topaz. The wooden surround is in beech finish/solid beech. Hand control with selective locking function.

manufactured according to:

- > norm DIN EN 1970:2000 + A 1:2005
- > norm DIN EN 60601-2-38/A1:2001 (reduced for home care requirements)
- > norm IEC 60601-2-52:2009

> order numbers		> order numbers		> dimensions
51.3431.33\$	home care bed "Arminia III" (wooden side rails)	51.3424.33	home care bed "Arminia III", 24-volt drive system (wooden side rails)	length 220 cm, width 101 cm
19.1400.1	lifting pole (included)	19.1400.1	lifting pole (included)	

> supplementary equipment (available at extra cost):

(technical changes reserved)



serving tray



mattress base with metal slats



clamp-on extra-high side rail



reading lamps




bed extension with wooden slats or metal slats (4 longer wooden side rails 220 cm) (also for retrofitting)



padded side rail cover

Arminia III - the practical home care bed in a homely design

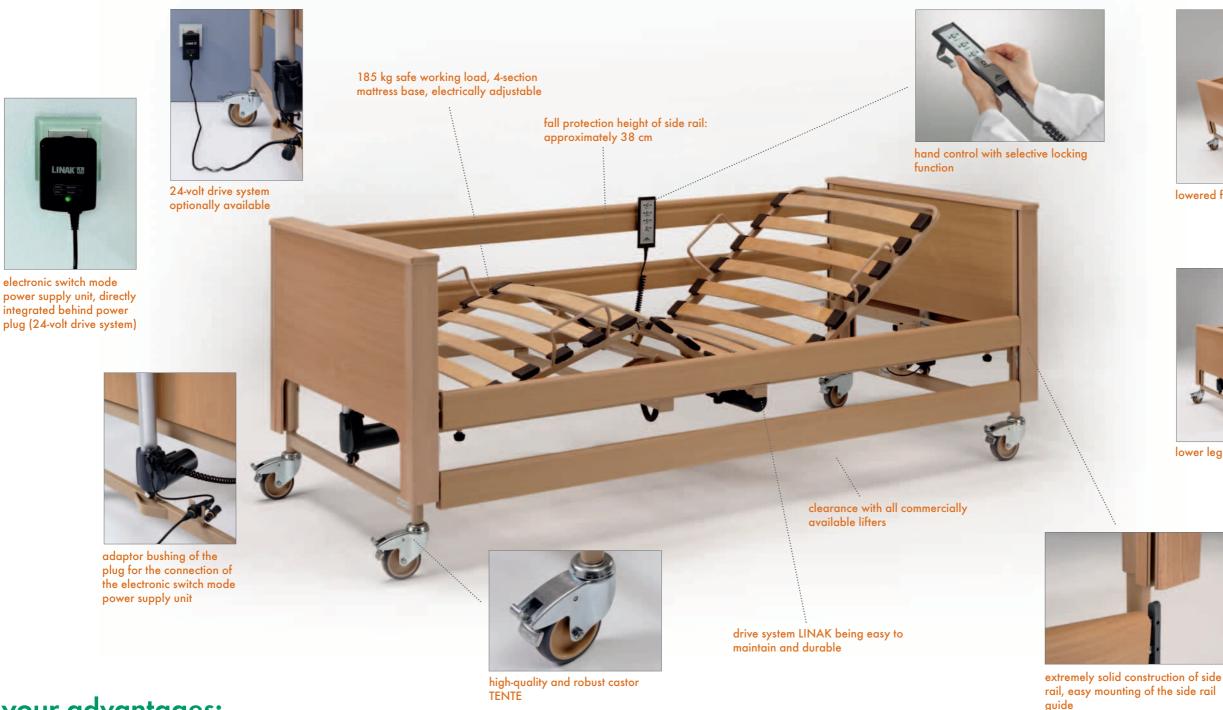



quality and durability

The home care bed Arminia III is extremely solid constructed, so that it fulfills the requirements in daily care perfectly. The high quality ensures you a home care bed which is durable, robust and easy to maintain and which decisively minimizes the follow-up costs. The wooden surround creates a homely ambience. A bed in which the resident feels very comfortable.

economic efficiency

The price-performance ratio is exceptionally reasonable. Nevertheless there was no economizing concerning the construction and the material. Especially the use of highquality materials ensures Burmeier's quality philosophy.

• 24-volt drive system optionally available

Also with this system the maintenance intervals can be extended to 2 years (see user manual). Leakage current measurements at the home care bed in the interval of 10 years. Conditions for this can be enquired from Burmeier. Residents are protected in the best possible way, as all 230 volt powered components are outside of the bed.

• fall protection height of side rail: approx. 38 cm Mattress and pressure ulcer systems of max. 16 cm thickness can also be placed in the home care bed without adaptable side rail raisings. In many cases costly raisings are not necessary.



lowered foot section



lower leg rest liftable by ratchet