
# Dali II 24 Volt



Clinic | Care | Homecare | Assist



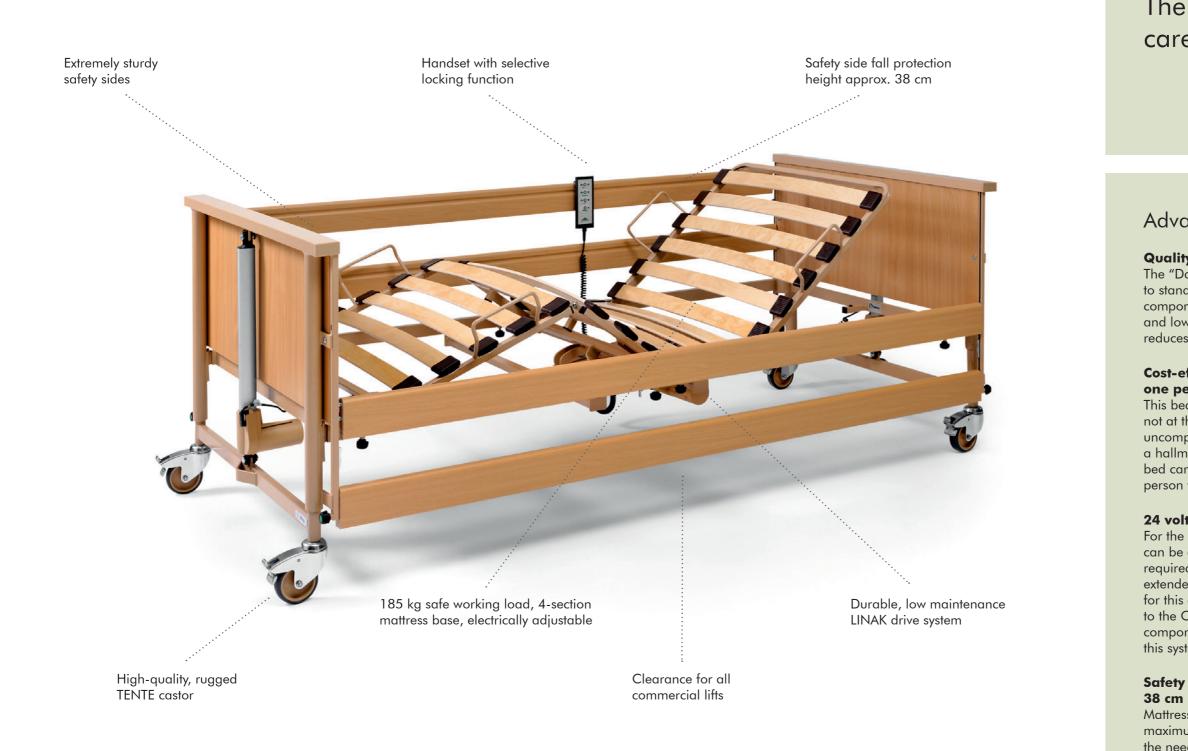

### Design details



Reverse-Trendelenburg position



Dali II with metal safety sides (model 15.0624.11)



Electronic switch mode power supply unit integrated immediately downstream of the mains plug



Power socket for connecting the electrical power supply unit (switch mode)

# **Dali II 24 Volt** -The cost-effective care bed

### Advantages at a glance

#### **Quality and durability**

The "Dali II 24 Volt" care bed is sturdily constructed to stand up to the rigours of daily care. High-quality components make this a particularly durable, sturdy and low-maintenance care bed that significantly reduces follow-on costs.

# Cost-effective / Easily assembled by just one person

This bed offers exceptional value for money, but not at the expense of materials or design. The uncompromising use of high-quality materials is a hallmark of the Burmeier philosophy. The care bed can be assembled and disassembled by one person without tools in a matter of minutes.

#### 24 volt drive system available as standard

For the 24 volt drive system, the maintenance interval can be extended to 2 years, and the interval for the required leakage current measurement can even be extended to a maximum of 10 years. The conditions for this can be obtained from Burmeier – please refer to the Operating Instructions. Since all live 230 volt components are located separately from the care bed, this system offers optimum safety for the occupant.

# Safety side fall protection height approx. 38 cm

Mattresses and anti-decubitus systems up to a maximum thickness of 16 cm can be used without the need for adjustable safety side height extensions. In many cases costly height extensions are not required at all.



Designed for space-saving transportation and storage

# Dali II 24 Volt care bed

### Model features

- Mattress base: 90 x 200 cm, 4-section, mattress base height adjustable from approx. 40 80 cm, backrest, thigh rest and bed height adjustable with electric motors
- Reverse-Trendelenburg position adjustable by tilting the mattress base
- Metal parts finished with polyester powder coating in the colour topaz
- Headboard and footboard panels in beech decor
- Handset with selective locking function
- Safe working load: 185 kg

Manufactured to: standard IEC 60601-2-52:2009

Accessories

(available at extra cost)



Mattress base with metal slats, manually adjustable lower leg rest



Wooden panels to cover the chassis (also for retrofit\*)



Bed extension including safety side bars (also for retrofit)



Reading lamps



Add-on safety side bar



Foam leather cover

### Order numbers and dimensions

- 51.0624.33 Care bed **Dali II 24 Volt** (wooden safety sides), length 218 cm, width 101 cm
- 15.0624.11 Care bed Dali II 24 Volt (metal safety sides)
- 19.1400.1 Patient lifting pole (included)

 Burmeier GmbH & Co. KG

 Pivitsheider Straße 270, 32791 Lage, Germany

 Phone:
 +49 (0) 5232 9841 - 0

 Fax:
 +49 (0) 5232 9841 - 41

 E-mail:
 info@burmeier.com

 www.burmeier.com
 www.burmeier.com



CE